




Maya

limbs&
things

Labour & Birthing Simulator

80500

80600

Fully automated birthing to teach the 3 stages of labour.



limbsandthings.com



Fully automated birthing to teach the 3 stages of labour.

Automated training for vaginal examination, delivery scenarios and post-partum haemorrhage. Set-up the model in minutes, select your training parameters and let Maya do the rest.



Quick to set-up and use so you can focus on what matters most



Repeatable training means consistency for every student



Quality teaching time: focus on teaching with fully automated simulation



Independent learning: ease of use means students can learn in their own time



Suitable for use with a standardised patient for added realism



Mess free post-partum haemorrhage



Ideal for examinations and OSCEs



No subscription required



Maya's 3 modes

VAGINAL EXAMINATION

Practice examination of the active stages of labour with realistic landmarks including ischial spines and vaginal wall.


CONFIGURATIONS

- 4 presentations:
 - Flexed
 - Deflexed
 - Indifferent
 - Breech
- 4-10 cm dilation
- Full 360° positioning



DELIVERY

- Build fully automated birthing scenarios to suit your training sessions.
- Cephalic & breech
- Log student actions as the scenario progresses & receive a time stamped report
- Sensors automatically detect McRoberts & suprapubic pressure interventions
- Can support forceps
- Clamp & cut umbilical cord
- IM injection site

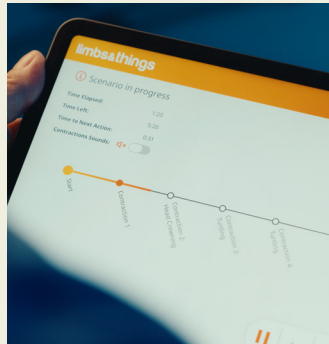

CONFIGURATIONS

- Foetal position:
 - OA
 - OP
 - Breech
- Select number & frequency of contractions
- Complications:
 - Shoulder dystocia
 - Assisted vaginal birth (ventouse, forceps & OdonAssist)
 - Post-partum haemorrhage

POST-PARTUM HAEMORRHAGE (PPH)

Available as a stand-alone training mode.

SCENARIO

Simulate full PPH in an integrated contained mat for a mess free bleed and easy clean up.


Configurations:

- Time onset after birth
- Blood flow rate
 - Mild - Moderate - Severe
- Intervention that will stop bleeding
 - Fundal massage
 - Removal of placenta piece
 - Manual stop

FUNDAL MASSAGE

As pressure is applied, the uterus will change between boggy (atonic) and contracted states. Perfect technique and get used to detecting the feel of each state.



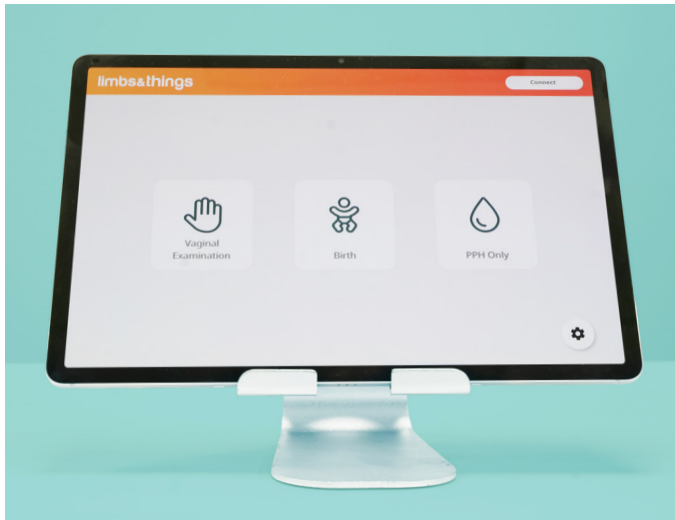

USER INTERFACE

- Maya is controlled using a simple app which connects to the model via Bluetooth & is available on any tablet.
- Save birthing scenarios in the app to make future set up consistent & quick.
- Advanced settings allow users to set preferences for position of the baby at certain milestones such as crowning.
- Available in 12 different languages.


MAYA'S BABY

- Representative of a full-term infant weight & size
- Spinal flexion
- Correct palpable landmarks (clavicles, scapula, suture lines & fontanelles)
- Fully articulating arms & legs

WHAT'S INCLUDED

- | | | |
|------------------------------------|-------------------------------|----------------------|
| • Maya | • 4 presenting parts & cervix | • Wheeled Case |
| • Baby | • PPH uterus | • Battery* |
| • Attachable lower legs | • Clean Bleed Mat & reservoir | • Charging cable |
| • Placenta (with retainable piece) | • Bed strap | *8 hour battery life |

Find out
more on our
website

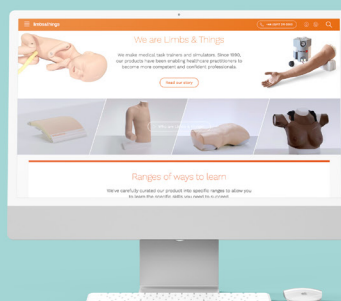

REQUEST A QUOTE



REQUEST A DEMO



BUY ONLINE



CUSTOMER SUPPORT

AUS V.1 August 2026

limbsandthings.com